

CELLNERGY WELLNESS



INSTRUCTION MANUAL

[ENGLISH](#)

CLICK BELOW FOR OTHER LANGUAGES

[ČEŠTINA](#) | [DEUTSCH](#) | [MAGYAR](#) | [POLSKI](#) | [ROMÂNĂ](#) | [ESPAÑA](#) | [РУССКИЙ](#)



PRODUCT SPECIFICATION

Voltage: 110-127V for North America

220-240V for EU

Frequency: 50/60Hz

Power: 900W

Sound: 80dB

SAFETY INSTRUCTIONS

Please read the instruction manual carefully before using.


- Keep device away from water.
- Please supervise children under 12 years old or anyone with physical or mental impairments who are using, cleaning, and maintaining this device. Do not leave children unattended with device or allow children to play with device.
- Do not use device near baths, showers, sinks, or near any liquids.
- Please be sure to unplug the device after turning it off after each use, especially in bathrooms.

WARNING: ELECTROCUTION HAZARD

Do not submerge this device in water or other liquids. If the device comes into contact with water, immediately disconnect it from any power source and do not touch it. Contact a qualified electrician for assistance. Failure to follow these instructions may result in electrocution or other serious injury. Keep this device away from children and water sources and use only in a dry environment. Do not use this device if it has been damaged by water or other liquids.



DANGER: KEEP AWAY FROM WATER!

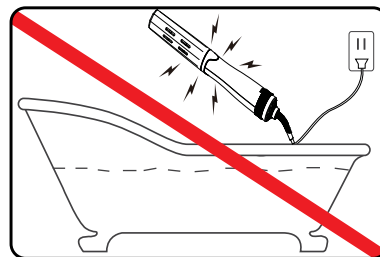

CAUTION

1. Read instruction manual carefully before use.
2. Do not allow children to play with device.
3. Please avoid knocking or dropping device.
4. Do not use the device outdoors, it should only be used in households.
5. Remove contact lenses before using device on face.
6. Do not plug in device if setting is not switched to OFF beforehand.
7. Do not blow device for more than 3 seconds on skin without movement.
8. Do not block the airflow of the device.
9. Do not use the device for more than 15 minutes per treated area daily.
10. Turn the setting to OFF before unplugging the device after each use.
11. Do not immerse the cord, plugs, or device in any liquids.
12. Do not place or store device where it can be fallen or drop into a tub or sink.
13. Do not reach for an appliance that has fallen into water. Unplug it immediately.
14. Do not use while bathing or in the shower.
15. Do not use the device if the supply cord is damaged.
16. If the supply cord arrived damaged, please contact our customer service team.
17. Do not leave the device unattended.
18. Unplug this appliance before cleaning.

SAFETY GUIDE

- Keep the product out of reach of children.
- Close supervision is necessary when the appliance is used by, on, or near individuals with certain disabilities or special needs.
- Do not leave appliance unattended when plugged in.
- Do not disassemble or modify the product in any way.
- Do not use the product near water or in wet conditions.
- Do not place the product near an open flame or any sources of heat.
- Do not expose the product to extreme temperatures or direct sunlight.
- Make sure the voltage of the product corresponds to the voltage of your power outlet before plugging in the device.
- Only use the product with the supplied power adapter.
- Do not use the product if the power adapter or cable is damaged, if it has been dropped or damaged, or if it has been submerged in water.
- Disconnect the device by gently pulling the plug of the power cord, not by pulling on the power cord directly.
- Do not use the product while driving or operating heavy machinery.
- Do not use the product if you have a pacemaker or any other medical device that could be affected by electromagnetic fields.
- Disconnect the device from the power socket after each use.
- Do not repair or make any adjustments to the device.
- Keep cord away from heated surfaces. Do not wrap the cord around the appliance.
- Do not use while sleeping.
- Do not operate when aerosol (spray) products are being used, or where oxygen is being administered.
- Do not use extension cord.
- Always store device in moisture-free area.
- Do not store device in temperatures exceeded 140°F (60°C).
- Do not unplug or plug this device when your hands are wet.

USER GUIDELINES

1. Do not use this product if you have hemorrhagic disease, a serious infection, hypertension, or open wounds.
2. If you have heart disease or any other medical condition, please check with your physician before using this product.
3. Do not use this product if you are undergoing surgical operations or pregnancy.

Disclaimer: This product is used for self-education and research purposes to support homeostasis. It is not intended to treat disease or to prevent impairment of human health. When used according to directions, it does not present any potential risk of illness or injury.

INSTRUCTIONS

1. Familiarize yourself with the settings:

Setting 0: Off

Setting 1: High speed, low heat

Setting 2: Low speed, high heat

Setting 3: Medium speed, high heat



We recommend Setting 1 for facial treatments, Setting 2 for ribcage and spine for stem cell activation, and Setting 3 for targeted areas.

2. Make sure the setting is OFF before plugging the device to the power outlet.
3. Plug the device into a power outlet and then use the dial at the bottom to adjust the setting.
4. Hold the device 8 to 12 inches (20 to 30 centimeter) away from the body. You can either blow the hot air in sweeping back and forth motions on a large area over bare skin or over clothing. Spend 10 minutes on each treated area on the body.
5. When you have finished using the device, turn it off before unplugging it from the power outlet.

WARRANTY

One-year warranty (from date of purchase) for product. If the device is defective during the warranty period, please contact our Customer Service team through our website www.lifepharm.com or by emailing customerservice@lifepharm.com for a new or refurbished replacement.

Warranty will be void if:

1. Warranty period is exceeded.
2. The device is damaged due to not following the instruction manual.
3. The device is damaged due to mishandling or repair by an unauthorized party.

60 DAY MONEY-BACK GUARANTEE


Returning product for a refund is easy! LifePharm offers a “no-hassle” return policy for retail customers. Committed to your complete satisfaction, we gladly back all our products with a 60-day money-back guarantee. All orders can be returned within 60-days from the order delivery date. All returns must be postmarked within the 60-day guarantee period. The cost of a product(s) is refundable, while shipping and handling are not refundable. We charge a \$30 handling fee for returning Cellnergy Wellness, and return postage is the buyer’s responsibility. Please visit our website or contact our customer service team for more information on our return policy.

MAINTENANCE

1. Always unplug and allow device to cool before cleaning.
2. Periodically inspect the air filter screen on the back of the device and remove any debris that might have accumulated. Any build up may reduce air flow and damage the device.
3. Never use thinner, gasoline alcohols, or any flammable chemicals to clean the product.
4. Clean the device with a soft dry cloth to remove any residue.
5. Do not wrap the cord around the device, this might cause the cord to wear out.

STEPS OF USING CELLNERGY WELLNESS

- Drink water before and after using the device.
- Hold the device 8 to 12 inches (20 to 30 centimeter) away from the body.
- Blow the hot air in sweeping back and forth or circular motions over bare skin or over clothing. Spend 10 minutes on each treated area on the body.

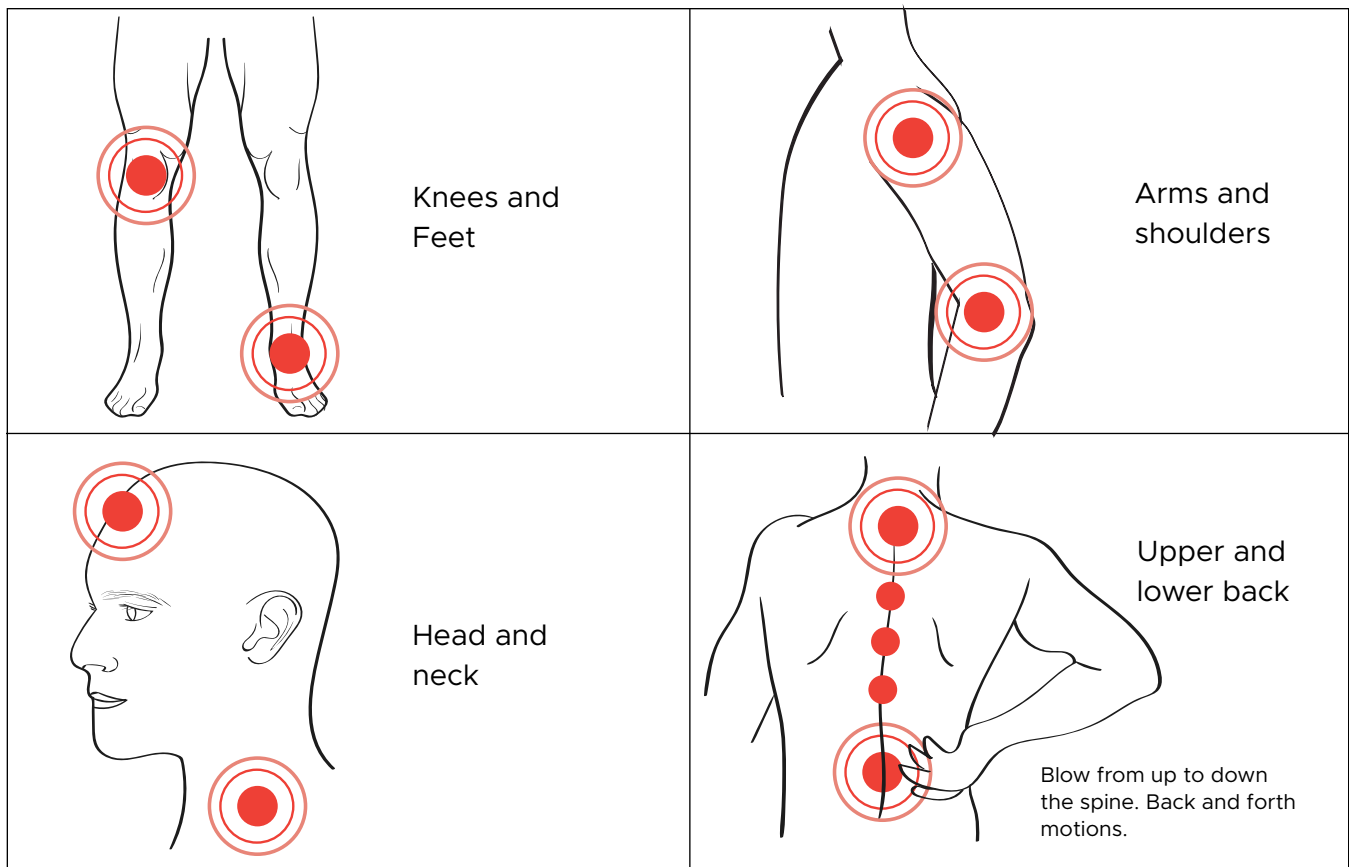

We recommend Setting 1 for facial treatments, Setting 2 for ribcage and spine for stem cell activation, and Setting 3 for targeted areas.

For best results, we recommend using Cellnergy Wellness daily with LAMININE and OMEGA+++ supplements.

FAQs

Q: What settings does the device have?

- A: The device has 4 settings:
- Setting 0: Off
 - Setting 1: High speed, low heat
 - Setting 2: Low speed, high heat
 - Setting 3: Medium speed, high heat

Q: What is the warranty?

- A: We offer a one-year warranty from date of purchase. If your device is defective during the warranty period, please contact our Customer Service team.

The warranty will be voided if the warranty period is exceeded, the device is damaged due to not following the instruction manual, or if the device is damaged due to mishandling or repair by an unauthorized party.

Q: Where is the product made?

- A: Our device is manufactured in China.

Q: Does the device have any kind of crystal inside?

- A: Yes, the device has a quartz crystal inside so the red light can pass through it safely due to its optical transparency, birefringent properties, piezoelectric effect, UV filtering capabilities, high-temperature stability, and low thermal expansion.

Q: What is the source of terahertz waves in the Cellnergy Wellness device?

- A: The terahertz waves in the Cellnergy Wellness device originate from a terahertz chip integrated into the device.

Q: Can terahertz waves be seen by the human eye?

- A: No, terahertz waves cannot be detected by the human eye because they are much longer in wavelength than visible light waves. They do not emit any visible light color.

Q: What is the voltage of the product?

- A: For North America, the voltage of our device is 110-127 V. For Europe, the voltage of our device is 220-240 V.

Q: Does the device restructure water?

- A: Terahertz waves have been proven to restructure water molecules to a hexagonal shape. If you would like to restructure water, we recommend waving the device over your glass of water for 30 seconds to 2 minutes.

Q: Will the device overheat if left on for too long?

- A: We do not recommend using the wand for more than 10-minute intervals. After each 10-minute interval is over, we recommend turning the wand off for at least 5 minutes to cool off.

Q: Can I use the device if I have a preexisting health concern or condition?

- A: If you have any preexisting health concerns or conditions, we recommend consulting with your health provider first before using the device.

Q: Can I use the device if I am pregnant or menstruating?

- A: No, we do not recommend using the device if you are pregnant or menstruating, since it helps increase blood flow.

Q: Can I use the device on my scalp?

- A: Yes, the device can be used 8-12 inches away from the scalp if there are no open wounds present.

Q: Can the device be used on any metal or implants in the body?

- A: No, we do not recommend using the device near any metal or implants on the body.

Q: Can the device be used on any open wounds?

- A: No, the device should not be used on any open wounds on the body.

Q: If I have a skin condition like psoriasis or eczema but with no open wounds on the skin, can I use the device?

- A: Yes, the device can be used on skin conditions if no open wounds are present.

Q: Can I use the device when I am wearing electronics?

- A: We recommend taking off all electronic devices before using the device so they do not overheat.

Q: Does the device emit any Electromagnetic (EMF) Radiation?

- A: No, the device has been tested with an Electromagnetic Radiation Detector to not emit any radiation.